


Bead Dispensing Head for Nanodrop Stages



Innovadyne introduces the **Bead Dispensing Head and Harness** for reliable dispensing of SPA and other beads.

The Bead Dispensing Head converts any Nanodrop stage to a bead dispenser in a matter of minutes, and is just as easily removed to convert back to a conventional Nanodrop stage.

When used in conjunction with bead dispensing protocols in the Nanobuilder software, Innovadyne's Nanodrop platforms equipped with the bead-dispensing-head can deliver unsurpassed bead dispensing performance, plate after plate, even with difficult beads such as Yttrium Oxide.

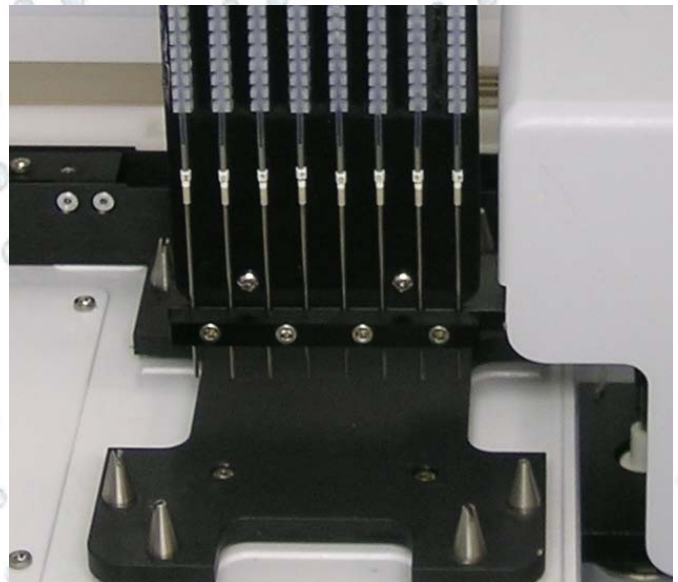

Features

- Bead suspensions are aspirated into vertical tubing, to minimize settling in the tubing and make it more uniform from tip to tip
- Quick installation and removal with hand tools, enabling the Nanodrop to be used with the conventional head for other applications and the bead dispensing head for bead dispensing
- Dedicated bead harness eliminates need to rethread harness into the head when changing heads

Ordering

| Item | Description |
|-------|--|
| 12087 | Bead dispensing head/harness assembly for Nanodrop I |
| - | Bead dispensing head/harness assembly for Nanodrop II/ExtY |